QUICK USER GUIDE

CHARGING THE BATTERY

The LifeinaBox battery needs to be fully charged before you start using the device

- 1) Make sure the device is connected to the power ①
- 2) Flick the O/I switch **under the device** to the ON position ②
- 3) Press the button **on the side** of the device ③ for 5 seconds until the logo in the front switches off.
- 4) Leave it to charge for two hours
- 5) Switch the device on by pressing the button on the side of device ③ and check that it is fully charged using the Lifeina application

USING LIFEINABOX


- Connect it to the power and fully charge the battery (see above)
- Switch on LifeinaBox by either:
 - 1) Pressing the button **on the side** of device once ③
 - 2) Unplugging the device and plugging it back in
- Wait for LifeinaBox to reach a stable temperature (about 20 minutes)
- Place your medication inside
- To use the Lifeina application, simply download the app and turn it on within distance of your LifeinaBox
- To switch off the logo on front of device, simply click the button on the side twice rapidly ③

NOTES

- If you are planning on using LifeinaBox at home, switch off the battery (under the device ②) so that it lasts longer. But don't forget to fully charge it before a trip
- If you want to use the battery mode, flick the O/I switch under the device to the "ON" position

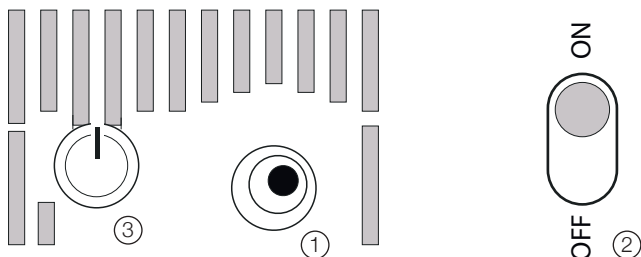Please do not block the air vents on the side and under the device

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